

UM-SJTU JOINT INSTITUTE

PHYSICS LABORATORY

VP141

Exercise 2

Measurements of Fluid Viscosity

Name: Fubo Qi
ID: 516370910177
Group: 11

July 4, 2017

I Introduction

I.1 objective

The objective of the exercise is to measure the fluid viscosity, which is one of the most important properties of fluids, determining the fluids flow, by using a common and simple method called Stokes' method.

I.2 Theoretical Background

When an object moves in a fluid, its motion is hindered by a drag force acting in the opposite direction of the movement of object. Also, the magnitude of the drag force is related to several components such as, the shape, speed of the object as well as the internal friction in the fluid. Usually, we use coefficient η to quantify the internal friction in the fluid, and the η is known as the viscosity coefficient. The drag force of a spherical object with radius R moving at speed v in an infinite volume of a liquid is given by,

$$F_1 = 6\pi\eta v R \quad (1)$$

And the buoyancy force is,

$$F_3 = \frac{4}{3}\pi R^3 \rho_1 g$$

Where ρ is the density of fluid and g is the acceleration due to gravity. The weight of the object is

$$F_3 = \frac{4}{3}\pi R^3 \rho_2 g$$

After sometime, the ball will achieve equilibrium speed when the three forces balance each other,

$$F_1 + F_2 = F_3. \quad (2)$$

Together with the equation 1, we get,

$$\eta = \frac{2}{9}gR^2 \frac{\rho_2 - \rho_1}{v_t}. \quad (3)$$

Because of the constant velocity, the η can also be written as

$$\eta = \frac{2}{9}gR^2 \frac{(\rho_2 - \rho_1)t}{s}, \quad (4)$$

where s is the distance traveled in time t with reaching the terminal speed.

Since the volume of the fluid is not infinite, the results are affected by some boundary effect

due to the presence of the container. Therefore the formula for the corrected magnitude is given by,

$$F_1 = 6\pi\eta v R \left(1 + 2.4 \frac{R}{R_c}\right)$$

Where R_c is the radius of the infinitely long cylindrical container. Consequently, we get,

$$\eta = \frac{2}{9} g R^2 \frac{(\rho_2 - \rho_1) t}{s} \frac{1}{1 + 2.4 \frac{R}{R_c}}. \quad (5)$$

II Apparatus and Measurement procedure

II.1 Experimental Setup

This exercise requires a Stokes' viscosity measurement device with castor oil and some small metal balls. There are also a number of measurement devices including micrometer, calliper, densimeter, electronic scales, stopwatch, and thermometer.

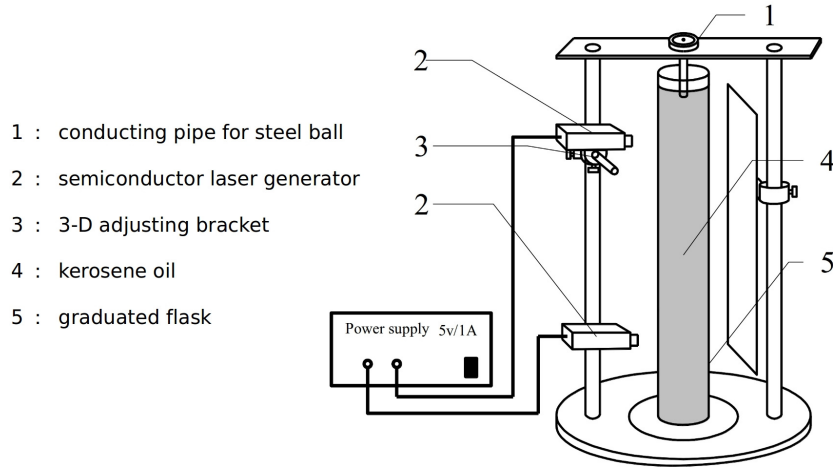


Figure 1: Stokes viscosity measurement apparatus

II.2 Measurement

II.2.1 Adjustment of the Stokes' viscosity measurement device

1. Make the plumb aiming at the center of the base by adjusting the knobs beneath the base.

2. Turn on the two lasers and make the beams parallel and aim at the plumb line.
3. Remove the plumb and place the flask with castor oil at the center of the base.
4. Place the guiding pipe on the top of the device.
5. Put a metal ball into it and check whether the ball can block the beams. If not, repeat.

II.2.2 Measurement of the (constant) velocity of a falling ball

1. Measure the vertical distance s between the two laser beams at least three times.
2. Put a metal ball into the guiding pipe. Record the time it travels between the two beams for at least six times.

II.2.3 Measurement of the ball density ρ_2

1. Use electronic scales to measure the mass of 40 metal balls. Calculate the average.
2. Use a micrometer to measure the diameter of the ball for ten times.
3. Calculate the density ρ_2 .

II.2.4 Other measurements

Measure of the density ρ_1 of the castor oil by using the provided densimeter (one measurement). Use a calliper to measure the inner diameter D of the graduated flask for six times. Read the ambient temperature from the thermometer placed in the lab.

III Results

III.1 Measurement of the Distance

<i>Distance</i> $S [\times 10^{-3}m] \pm 0.10 [\times 10^{-3}m]$			
S_1	110.0	263.0	153
S_2	140.0	292.0	152
S_3	111.0	264.5	153.5

Table 1: Distance measurement data

The average value of the vertical distance of two beams is hence calculated by,

$$\bar{S} = \frac{1}{3} \sum_{i=1}^3 S_i = (152.8 \pm) \times 10^{-3}m$$

title

III.2 Measurement of the Time

Measurement <i>Time</i> t [s] \pm 0.01 [s]	
t_1	6.97
t_2	7.12
t_3	7.28
t_4	7.35
t_5	6.94
t_6	7.00

Table 2: Time measurement data

The average value of the time is given as

$$\bar{t} = \frac{1}{6} \sum_{i=1}^6 t_i = 7.11 \pm s$$

III.3 Diameter of ball

<i>Diameters</i> $d[\times 10^{-3}m] \pm 0.004[\times 10^{-3}m]$			
d_1	1.945	d_6	1.955
d_2	1.945	d_7	1.940
d_3	1.950	d_8	1.940
d_4	1.955	d_9	1.950
d_5	1.950	d_{10}	1.935

Table 3: Measurement data for diameters of the balls

With the initial reading, $d_0 = -0.040mm$ The average value of the diameter of a ball is hence calculated as

$$\bar{d} = \frac{1}{10} \sum_{i=1}^{10} d_i = (1.986 \pm) \times 10^{-3}$$

III.4 Inner Diameter of the Flask

The average value of the inner diameter of the flask is hence calculated

$$\bar{D} = \frac{1}{10} \sum_{i=1}^{10} D_i = (61.69 \pm 0.21) \times 10^{-3} m$$

<i>Diameters $D[\times 10^{-3}m] \pm 0.02[\times 10^{-3}m]$</i>	
D_1	61.70
D_2	61.78
D_3	61.68
D_4	61.70
D_5	61.64
D_6	61.68

Table 4: Measurement data for the inner diameter of the flask

III.5 Other physicals quantifies

Density of the castor oil	$\rho_1 = 0.956 \pm 0.001 g/cm^3$
Mass of 40 metal balls	$m = 1.359 \pm 0.001 \times 10^{-3} kg$
Temperature in the lab	$T = 26 \pm 2^\circ C$
Acceleration due to gravity in the lab	$g = 9.794 m/s^2$

Table 5: Values of other physical quantities

IV Uncertainty Analysis